# Incubator

IN260

The incubator I is at home everywhere in the world of research, medicine, pharmaceutics and food analytics, as well as food chemistry.




**Experts in Thermostatics** 

E

The heating of this incubator is optimally tuned for both natural convection and forced air circulation; the fan can also be switched off completely, and valuable chamber loads for research, pharmaceutics, medicine and food chemistry are warmed up very carefully.

On this page, you can find all the essential technical data on our incubator. Our customer relations team will be pleased to help if you want further information. If you should require a customised special solution, please contact our technical specialists at <u>myAtmoSAFE@memmert.com</u>.

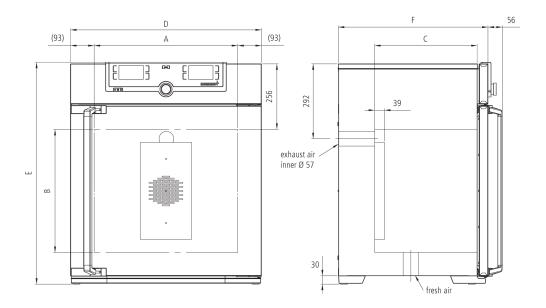

# **Control of standard components**

ControlCOCKPIT	adaptive multifunctional digital PID-microprocessor controller with high-definition TFT-colour display
Temperature	1 Pt100 sensor DIN class A in 4-wire-circuit
Timer	Digital backwards counter with target time setting, adjustable from 1 minute to 99 days

# Temperature

_	Set temperature range in °C	min. 5°C above ambient up to +80°C
	resolution of display for setpoint and actual temperature values	0.1°C

# **Control technology**

adjustable parameters	temperature (Celsius or Fahrenheit), air flap position, programme time, time zones, summertime/wintertime
Function SetpointWAIT	the process time does not start until the set temperature is reached
Language setting	German/English/Spanish/French
Calibration	three freely selectable temperature values

### Ventilation

natural convection	
adjustment of pre-heated fresh air admixture by air flap control in 10 % steps	
vent connection with restrictor flap	

Programming	AtmoCONTROL software for reading out, managing and organising the data logger via Ethernet interface (temporary trial version can be downloaded). USB stick with AtmoCONTROL software available as accessory (on demand).
Documentation	programme stored in case of power failure

# Safety

Temperature control	adjustable electronic overtemperature monitor and mechanical temperature limiter TB, protection class 1 according to DIN 12880 to switch off the heating approx. 20°C above nominal temperature
Autodiagnostic system	for fault analysis

Door	fully insulated stainless steel door with2-point locking (compression door lock)
Door	inner glass door
Internals	2 stainless steel grids
Scope of delivery	incl. works calibration certificate for +37°C
Housing	rear zinc-plated steel
Interior	easy-to-clean interior,made of stainless steel,reinforced by deep drawn ribbing with integrated and protected large-area heating on four sides

# Standard equipment

# **Stainless steel interior**

w <sub>(A)</sub> x h <sub>(B)</sub> x d <sub>(C)</sub> : 640 x 800 x 500 mm	
Volume 256 l	
Max. loading of chamber: 300 kg	

#### Textured stainless steel casing

 $w_{(D)} \ge h_{(E)} \ge d_{(F)} \ge 824 \ge 1183 \ge 684$  mm

#### **Electrical data**

Voltage	230 V, 50/60 Hz
Electrical load	approx. 1700 W
Voltage	115 V, 50/60 Hz
Electrical load	approx. 900 W

### Packing/shipping data

the appliances must be transported upright	
Customs tariff number	8419 8998
Country of origin	Federal Republic of Germany
WEEE-RegNo.	DE 66812464
	Dimensions approx incl. carton B x H x T: 930 x 1380 x 930 mm
	Net weight approx. 110 kg
	Gross weight carton approx. 161 kg

#### Standard units are safety-approved and bear the test marks

